

MPA-10/MPA-20/MPA-200/MPA-1200

Single Channel Electronic Pipette




INSTRUCTION MANUAL



1WMPD4002874B

WARNING DEFINITIONS

The warnings described in this manual have the following meanings:

 DANGER	An imminently hazardous situation which, if not avoided, will result in death or serious injury.
 WARNING	A potentially hazardous situation which, if not avoided, could result in death or serious injury.
 CAUTION	A potentially hazardous situation which, if not avoided, may result in minor or moderate injury or damage to the instrument.
NOTE	Information or cautions to use the device correctly.

© 2014 A&D Company Ltd. All rights reserved.

No part of this publication may be reproduced, transmitted, transcribed, or translated into any language in any form by any means without the written permission of A&D Company Ltd.

The contents of this manual and the specifications of the instrument covered by this manual are subject to change for improvement without notice.

CONTENTS

1. FOR SAFE USE	3
1-1 Precautions on using the pipette	3
1-2 Precautions on handling the battery	3
2. INTRODUCTION	4
3. FEATURE	4
4. COMPLIANCE	5
5. MPA FUNCTION	7
6. PACKING CONTENTS AND NAME OF ITEMS	8
7. PREPARATION BEFORE USE	10
7-1 Installing the battery	10
7-2 Recharging the battery	11
7-3 Exchange selectable power plug	12
7-4 Before operating the pipette	12
7-5 Parts names	14
8. NAME AND FUNCTIONS OF DISPLAY AND KEYS	15
8-1 Display and functions.....	15
8-2 Key switches and functions.....	16
9. FUNCTION AND HOW TO USE	17
9-1 Standard mode (AUTO).....	17
9-2 Multiple dispensing mode (MD)	19
9-3 Mixing mode (MIX)	21
9-4 System setting mode (SYS)	23
9-5 Program setting mode	25
9-6 Reverse operation (Dispensing liquid that tends to remain in the tip)	25
9-7 Blowout function	27
9-8 Total discharge function	28
9-9 “Dispensing correction function” for multiple dispensing	28
10. PIPETTING FOR ACCURATE DISPENSING	28
11. CALIBRATING THE PIPETTE USING AN ELECTRONIC BALANCE	29
11-1 Volume calibration function (μL calibration function)	29
11-2 Resetting the volume calibration.....	30
11-3 Dispensing in a unit of weight (in mg unit).....	31
11-4 Weight calibration function (mg calibration function)	32
12. ADJUSTING HEIGHT OF THE TIP EJECTOR	34
13. STORAGE AND MAINTENANCE	34
13-1 Autoclave.....	35
14. TROUBLE SHOOTING	36
15. WHEN REQUESTNG REPAIRE	37
16. SPECIFICATION	38

17. LIST OF ITEMS SOLD SEPARATELY (DISPOSABLE ITEMS) 39
17-1 Stands and hanger.....39
17-2 Tips and accessories40
17-3 Disposable items (User replaceable)42
17-4 Others43

1. FOR SAFE USE

1-1 Precautions on using the pipette

DANGER

- This instrument is not an explosion proof instrument. Do not use the pipette in an environment where there is a risk of explosion, or use it for explosive chemicals that may cause explosion.
- When using potentially harmful solutions, such as infectious bacteria or viruses, radioactive substances that have a risk of exposure, or poisons, exercise extreme caution and follow all safety measures.

CAUTION

- When using a corrosive solution such as an organic solvent, confirm the chemical-proofness of the tip or pipette. When confirming material of the pipette, refer to “7-5 Parts names”.
- Do not attempt to disassemble or repair the pipette by yourself. Refer to “14. TROUBLE SHOOTING” when it appears that the pipette has a mechanical error.

1-2 Precautions on handling the battery

The MPA series use the high-density lithium-ion battery.

To prevent injuries or accidents due to a leaking battery, heat generation, fire or burst, and to ensure safe use, be sure to keep the manual on hand.

DANGER

- Do not dispose of the battery in fire, do not heat it, do not disassemble or modify it.
- Do not splash water on the battery, or do not keep the battery in a location at high temperature or high humidity.
- Do not allow battery contacts to contact metal. When keeping or carrying the battery, be sure not to allow the battery to contact metal.

WARNING

- Recharge the battery with pipette installed. The pipette can be used even when the battery is being recharged.
- When recharging is unsuccessful even after charging for the specified time (Five hours up to fully recharge), stop recharging the battery.
- Use only the battery supplied with the pipette. Do not use other batteries.

CAUTION

- Do not use a leaking battery.
- Because the battery body may become hot when using the pipette continuously for a long time, take care not to get burned when handling it.
- Should you get battery fluid from a leaking battery in your eye, immediately flush with copious amounts of water and seek immediate medical attention.
Should the liquid contact the clothes or skin, rinse immediately with copious amounts of water.

2. INTRODUCTION

Thank you for purchasing the MPA series electronic pipette. To ensure safe use of the product, be sure to read the manual thoroughly.

3. FEATURE

The MPA series is a high precision and performance electronic pipette that achieves operability without putting a burden on the hand.

This pipette is developed for the purpose to prevent RSI (Repetitive Strain Injury) which may occur when repeatedly using a manual pipette, and does not require any special skill so anyone can easily and accurately dispense the specified volume.


- Pipette is operated by merely pressing a key, the degree*1 of fatigue is 1/100 or less of when using pipette manually. (*1 Calculated by operating force and movement)
- It has an ergonomic design, fitting the hand for easy adjustments and operation.
- Using a lithium-ion battery enables usage for long periods of time.*2
(*2 Approximately 1,800 dispensing operations are possible on a full charge)
- Impact-absorbing pads adopted to fully protect against falling. (Patent has been applied for)

4. COMPLIANCE

Compliance with FCC Rules

Please note that this device generates, uses and can radiate radio frequency energy. This device has been tested and has been found to comply with the limits of a Class A computing device pursuant to Subpart J of Part 15 of FCC rules. These rules are designed to provide reasonable protection against interference when this device is operated in a commercial environment. If this unit is operated in a residential area, it may cause some interference and under these circumstances the user would be required to take, at his own expense, whatever measures are necessary to eliminate the interference. (FCC = Federal Communications Commission in the U.S.A.)

Compliance With EMC Directives of CE mark

 This device features radio interference suppression, safety regulation and restriction of Hazardous Substances in compliance with the following Council Directives

Council directive 2004/108/EC	EN61326	EMC directive
Council directive 2006/95/EC	EN61010-1	Low voltage directive
Council directive 2011/65/EU	EN50581	Restriction of the use of certain Hazardous Substances

The CE mark is an official mandatory European marking.

Please note that any electronic product must comply with local laws and regulations when sold or used anywhere outside Europe.



A & D Instruments Ltd. hereby declare that the following Weighing product conforms to the requirements of the council directives on ...

**Electromagnetic Compatibility (EMC) 2004/108/EC,
Low Voltage Equipment (LVD) 2006/95/EC amended by 93/68/EEC and
Restriction of the use of certain Hazardous Substances (RoHS) 2011/65/EU**

provided that they bear the CE mark of conformity.

Model/Series....MPA Series

Standards applicable:

EN61326-1:2013

Electrical equipment for measurement, control and laboratory use -EMC requirements Part 1: General requirements

EN-61010-1:2010

Safety requirements for electrical equipment for measurement, control, and laboratory use. General requirements

EN-61010-2-101:2002

Safety requirements for electrical equipment for measurement, control and laboratory use. Particular requirements for in vitro diagnostic (IVD) medical equipment

EN-50581:2012

Technical documentation for the assessment of electrical and electronic products with respect to the restriction of hazardous substances

CE Mark first applied 30th May 2014

Signed for A&D Instruments in Oxford England 18th June 2014

P. Argus
Managing Director

5. MPA FUNCTION

- The pipette has three modes where advantages of electromotion are utilized. (Refer to “9. FUNCTION AND HOW TO USE”)
 - Standard mode (AUTO)

This is for basic pipette operation. In this mode, the pipette aspirates once and then dispenses once.
 - Multiple dispensing mode (MD)

This is for dispensing liquid on a microplate, etc. In this mode, the pipette aspirates once and dispenses several times.
 - Mixing mode (MIX)

This is a useful operation when uniformly mixing different types of liquids. In this mode, the pipette repeats a cycle of aspirating and dispensing.
- User setting allows storage within the pipette of up to nine programs containing operating mode and dispensed volume. By reading them out when necessary, operation for setting again can be omitted. Settings from the prior use are stored in memory even with the power turned off.
- The pipette is equipped with the reverse operation suitable for dispensing a liquid that has a tendency to remain in the tip. (Refer to “9-6 Reverse operation (Dispensing liquid that tends to remain in the tip)”)
- The pipette also has “Dispensing correction function” (Patent applied for) with multiple dispensing to cancel errors due to backlash. It enables the dispensing of liquids precisely without difference due to operators. (Refer to “9-9 “Dispensing correction function” for multiple dispensing”)
- Various kinds of tips can be used. (The height of the tip ejector can be adjusted) (Refer to “12. ADJUSTING HEIGHT OF THE TIP EJECTOR”)
- Calibration (adjustment) of dispensed volumes is easy. (User CAL function). Even differences in dispensed amounts due to tip differences can be corrected. (Refer to “11-1 Volume calibration function (μ L calibration function)”) (Patent applied for)
- Dispensing by weight is also available. Refer to “11-3 Dispensing in a unit of weight (in mg unit)”) (Patent applied for)

6. PACKING CONTENTS AND NAME OF ITEMS

Confirm that the following contents are all included.

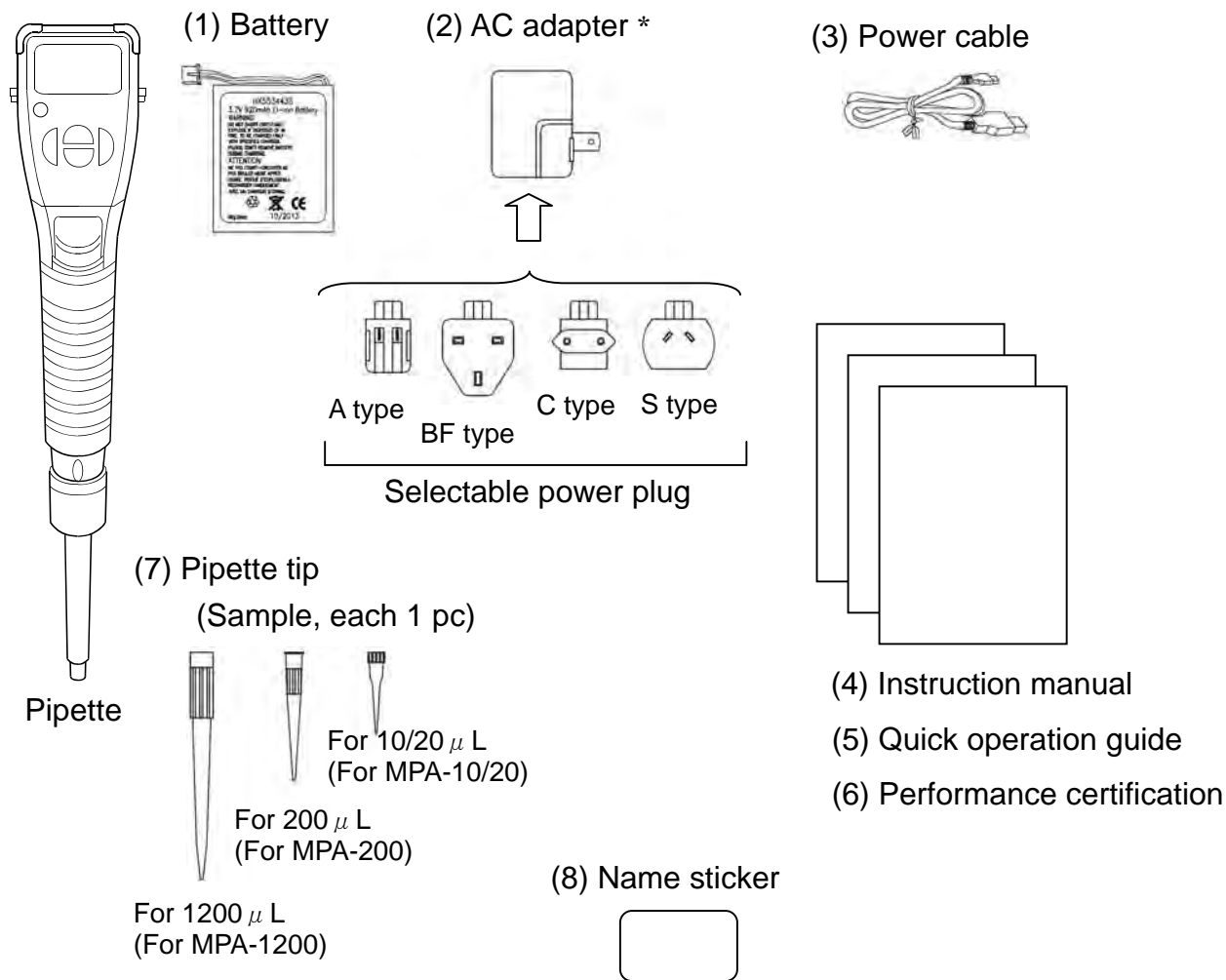
- Electronic pipette MPA-10 / 20 / 200 / 1200 (Any one among them)
- Accessories
 - (1) Battery (1 pc)
 - (2) The AC adapter (Combined use for charging) (Switching with AC100V to 240V)
Selectable power plug (A / BF / C / S type)
* AC adapter has the A type AC adapter plug attached.
Use AC adapter plug for AC adapter to match local outlet.

Note

Please confirm that the AC adapter type is correct for your local voltage and receptacle type.

- (3) Power cable (USB cable: Mini B plug - A plug)
- (4) Instruction manual (This document)
- (5) Quick operation guide
- (6) Performance certification (Pipette Accuracy Test Result)
- (7) Pipette tip (3 pcs)
(For 10/20 μ L (1 pc), for 200 μ L (1 pc), for 1200 μ L (1 pc))
- (8) Name sticker (The pipette has a location in the battery compartment for affixing the name sticker.

The pipette has a protective film on the display.
Remove it if necessary.



Should the pipette arrive damaged or an accessory be missing, contact the nearest A&D dealer.

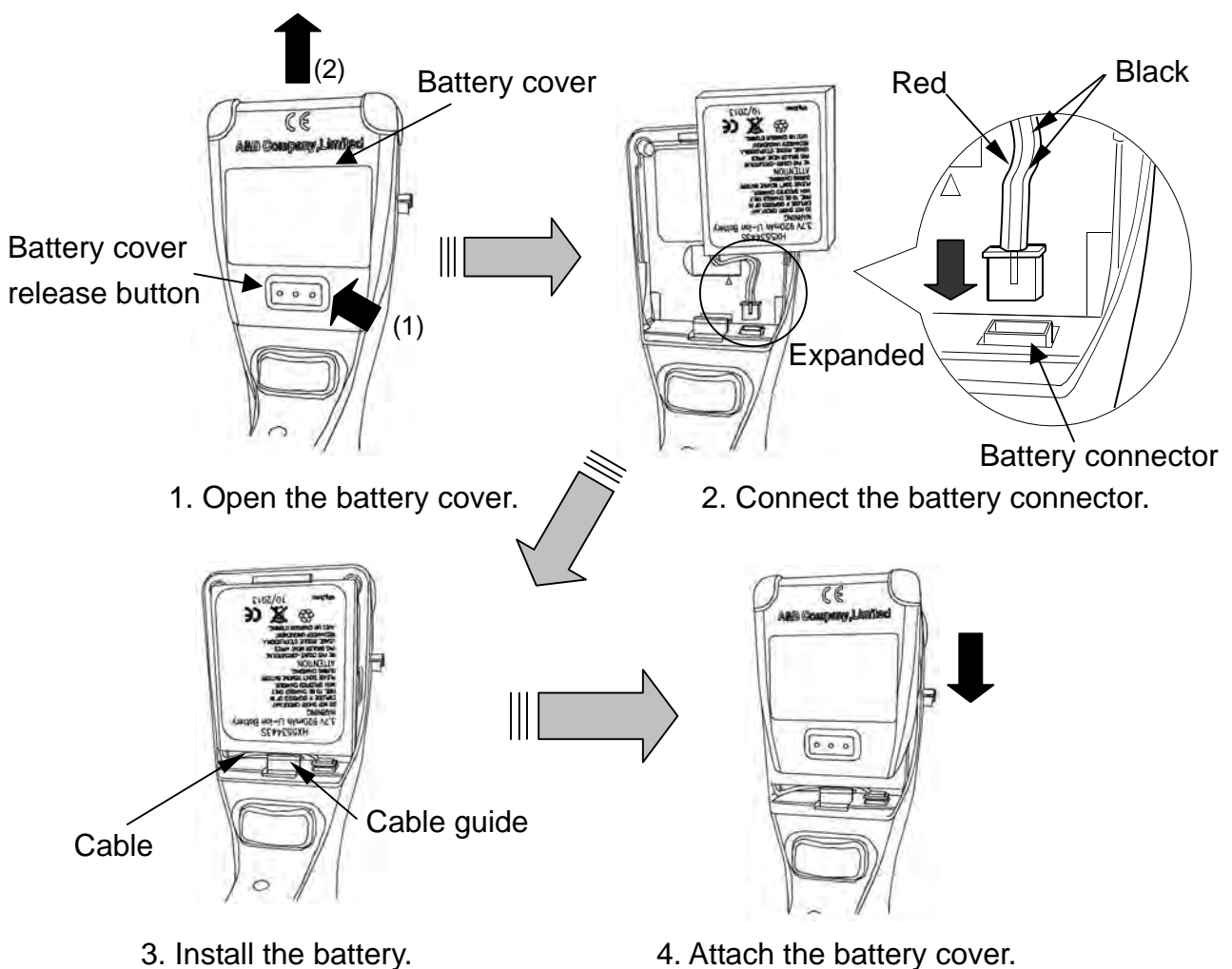
Note

The accessories included with this product may be changed without notice.

7. PREPARATION BEFORE USE

7-1 Installing the battery

1. Remove the battery cover (2) by sliding it upward while pressing and holding the battery cover release button (1).
2. Connect the terminal of the battery's cable, as shown in the figure below, to the connector for the battery in the bottom of the battery compartment. When connecting the terminal, be sure it is connected in the proper direction.
3. Install the battery so that the battery cable is in the cable guide.
4. Attach the battery cover on the pipette by sliding it downward from the upside.




Note

When connecting the battery to the pipette, all illuminations on the display illuminate and the pipette built-in piston automatically goes to the initial default position. If a key is pressed, the pipette goes into the operating mode.

7-2 Recharging the battery

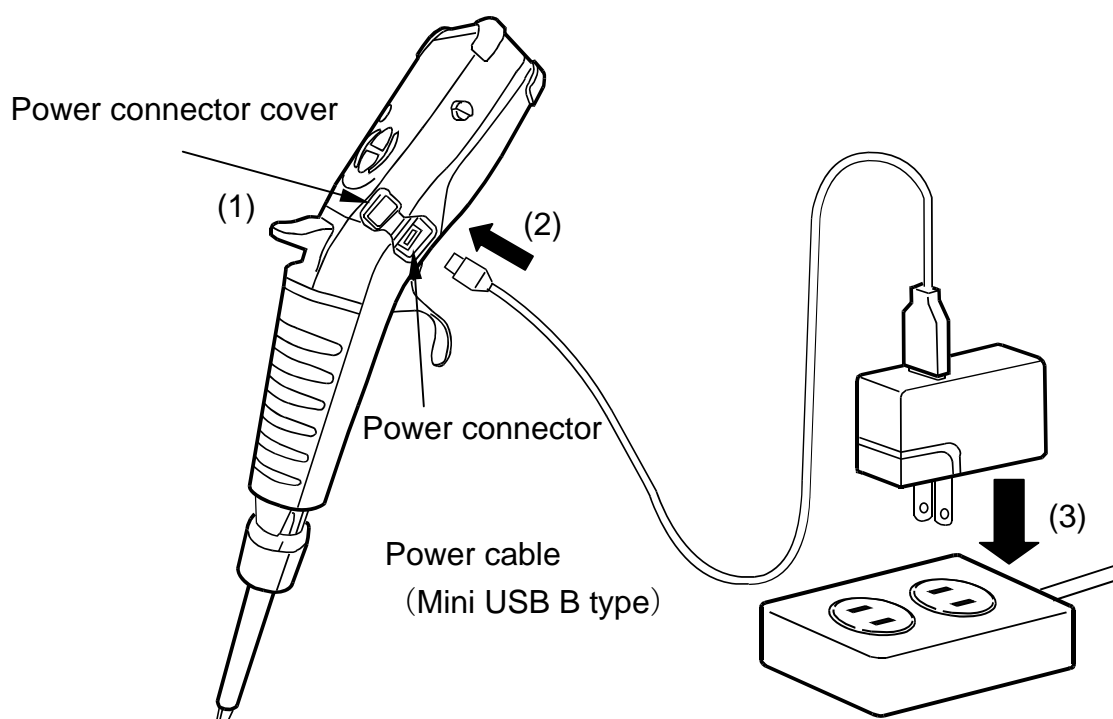
When purchasing the instrument, the battery does not have a full charge. For initial use, first charge the battery fully. Recharge the battery with the battery installed in the pipette. Pipette use is available during recharging.

Recharging

1. Remove the power connector cover from the pipette.
2. Connect the power cable connected to the AC adapter to the power connector on the pipette.
3. Connect the AC adapter plug to the outlet. The battery mark  will be displayed on LCD of the pipette, and it will blink during recharging. If connecting the power cable to the outlet before setting the battery in the pipette, please note that the recharging will not start. When the recharging is complete, the battery mark changes from blinking to a steady illumination, then the recharging completes automatically. (About five hours)

NOTE:

After recharging the battery completely, remove the power cable from the pipette. Firmly attach the power connector cover by pushing it onto the pipette.




7-3 Exchange selectable power plug

The A-Type power supply plug is originally attached to the AC adapter. Please change the power supply plug to the one that suits your location.

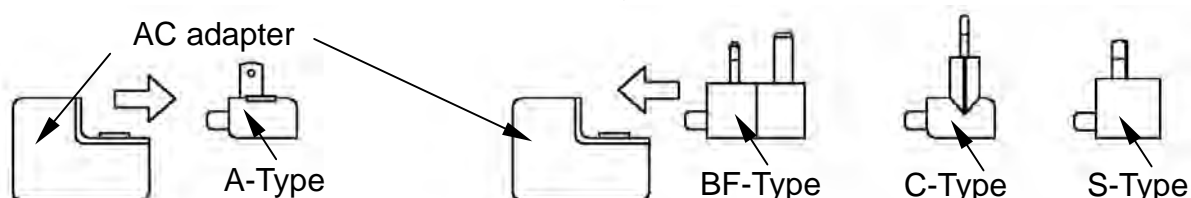Exchange method




1. As shown, remove the power supply plug from the AC adapter.
2. As shown, put on the power supply plug that you want to use.

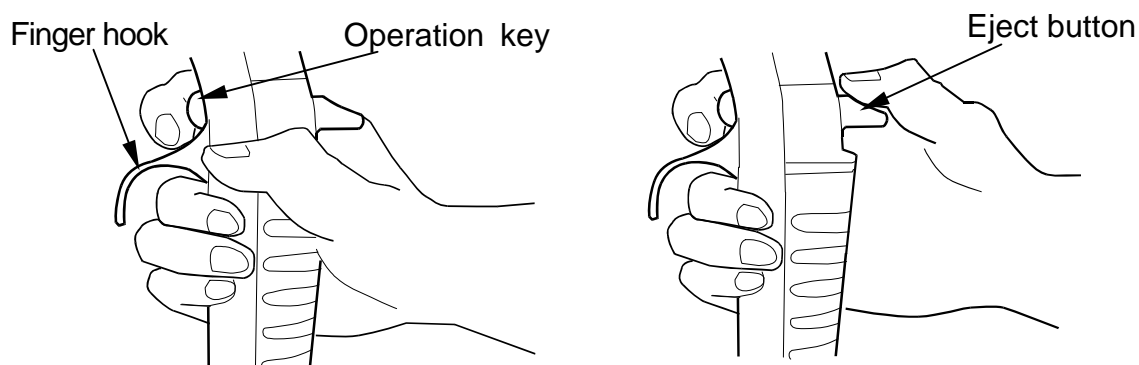


1. Remove the power supply plug.
2. Put on an appropriate power supply plug.


7-4 Before operating the pipette


Holding the pipette

- Hold the pipette so that the finger hook is between a forefinger and middle finger.
- To aspirate or dispense a liquid, operate the  Operation key or the  Up key below the display. Operate the  Operation key using the forefinger, as shown in the figure below.
- Operate the eject button by using the thumb to remove the tip.



Operating mode and standby mode

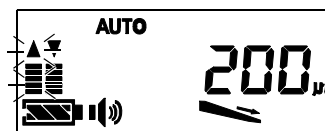
- The pipette goes into standby mode to reduce the battery wasting to minimize battery use if the pipette is idle for 10 minutes.
- When off, the pipette can be returned to the operation mode by pressing any key, and information such as setting volume will be displayed on the display (Refer to example of the display), enabling dispensing. At this time, the pipette automatically positions the built-in piston to the initial default position.
- While in the operating mode, holding down the  Operation key for approx. five seconds will turn the pipette off.





Example of display


Operation

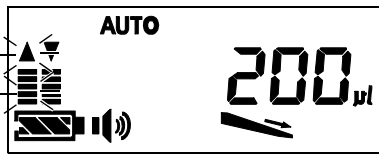Turning the power off manually




(The display is example.)

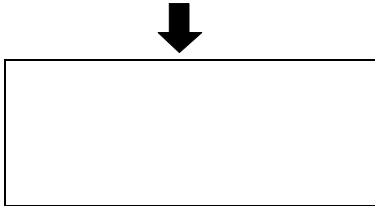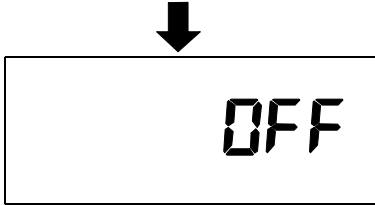

Press and hold the  Operation key on the pipette for at least five seconds until OFF is displayed on the display.

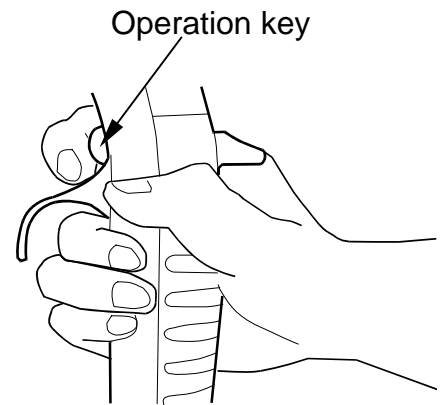
Buzzer sounds (Three times), and the pipette turns the power off (OFF).



7-5 Parts names


The following shows the each name of electronic pipette

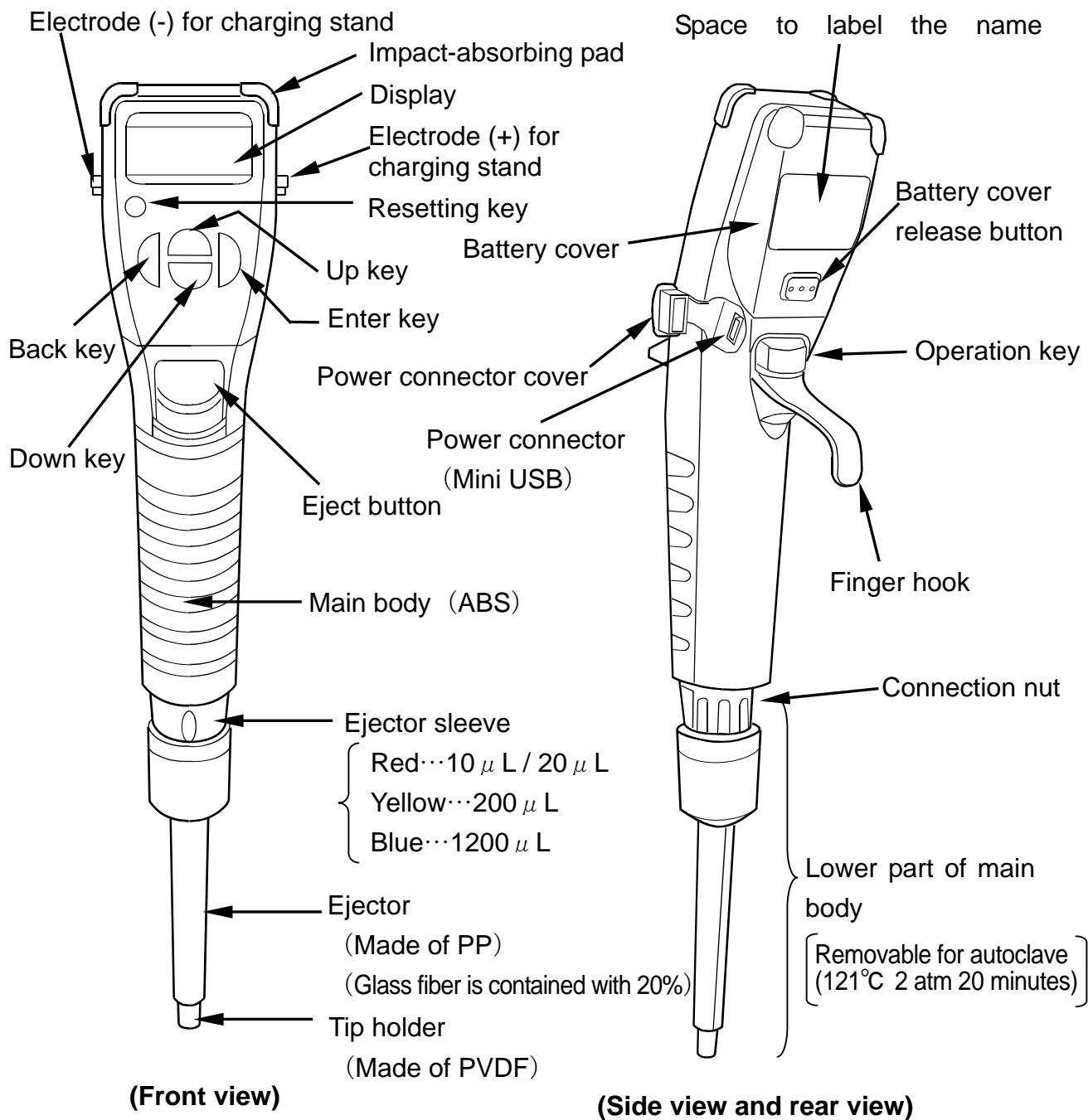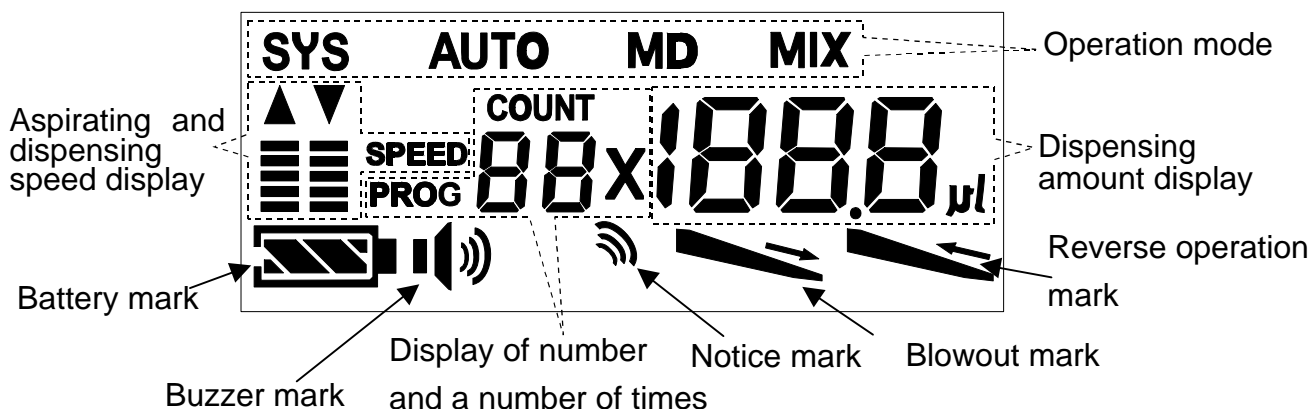When confirming LCD, refer to “8-1 Display and functions” for details.










8. NAME AND FUNCTIONS OF DISPLAY AND KEYS








8-1 Display and functions





	Symbols	Descriptions	
Operation mode display	SYS	System setting mode	Used to set up functions before pipetting, such as aspirating/dispensing speed, reverse operation, etc. (Refer to “9-4 System setting mode (SYS)”.)
	Shows the operation mode when operating the pipette.		
	AUTO	Standard mode	(Refer to “9-1 Standard mode (AUTO)”.)
	MD	Multiple dispensing mode	(Refer to “9-2 Multiple dispensing mode (MD)”.)
	MIX	Mixing mode	(Refer to “9-3 Mixing mode (MIX)”.)
Dispensing amount display	1888 µl	Shows the setting value of the dispensed amount. Shows the dispensed amount by unit of weight (Unit in mg) when not displaying the µL.	
Blowout mark		Shows whether the blowout is enabled or disabled. (Refer to “9-7 Blowout function”.)	
Reverse operation mark		Shows reverse operation. (Refer to “9-6 Reverse operation (Dispensing liquid that tends to remain in the tip)”.)	
Display of number and a number of time	COUNT 88x PROG 88x	COUNT : Shows the number of times the same operation is to be carried out. PROG : Shows the stored number of user setting. (Refer to “9-5 Program setting mode” .)	

	Symbols	Descriptions
Aspirating and dispensing speed display		Shows the speed level when aspirating or dispensing the liquid. ▲ blinks when aspirating, ▼ blinks when dispensing (Refer to “9-4 System setting mode (SYS)” .)
Notice mark		When illuminated: Shows that volume calibration has been carried out. When flashing: Shows that weight mode for dispensing (mg) has been selected. (Refer to “11-3 Dispensing in a unit of weight (in mg unit)”.)
Buzzer mark		Shows the buzzer is to sound or not. (Refer to “9-4 System setting mode (SYS)”.)
Battery mark		Shows the battery status.  Charging amount: Full  Charging amount: Low (Recharge the battery using AC adapter.)  During charging

8-2 Key switches and functions

Keys	Symbols	Functions and descriptions
Setting keys	Enter key 	Confirms the setting content.
	Back key 	Changes the mode or cancels it.
	Up key 	Increases the volume and setting value. Changes items (Mode).
	Down key 	Decreases the volume and setting value. Changes items (Mode).
Resetting key		Stops dispensing a liquid and returns the built-in piston to the initial default position. By pressing the Resetting key, all illuminations illuminate. After that, the pipette returns to the operating mode if any key is pressed.
Operation key		Starts aspirating and dispensing. Discharges all the liquid left in the tip when held down in the middle of multiple dispensing. Puts the pipette in standby mode when held down further.
Eject button		Detaches the tip.

Useful use method: The  Operation key (key switch on rear side on the pipette) has the same function as the  Up key. This allows you to quickly perform settings such as changing the volume without shifting the pipette in the hand.

9. FUNCTION AND HOW TO USE

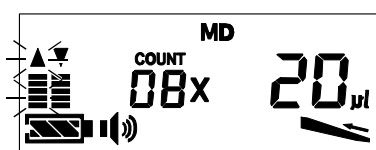
The MPA series have three modes, the standard mode (AUTO), multiple dispensing mode (MD) and the mixing mode (MIX).

9-1 Standard mode (AUTO)

1) Operating the standard mode


This is a basic operation for pipetting. Aspirating one time and dispensing one time.
This operation is the same as for a manual pipette.

2) Selecting the standard mode





(The display is example.)

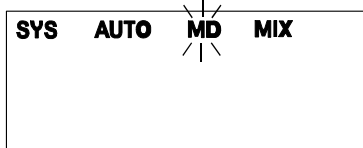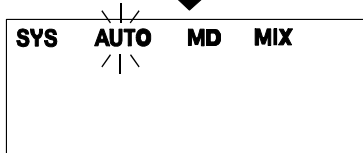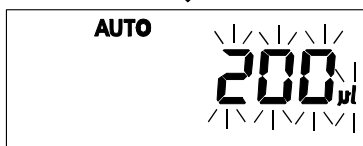
[1] Press the  Back key.




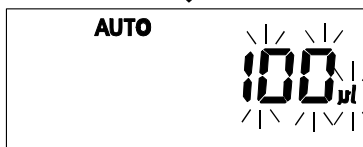
[2] Press the  Up or Down key to select "AUTO".




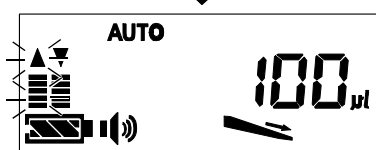
[3] Press the  Enter key to select standard mode (AUTO).




[4] Press the  Up or Down key to change to volume that you would like to dispense.



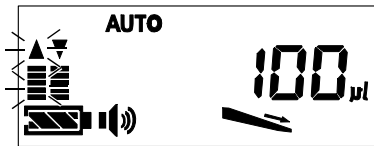
[5] Press the  Enter key to confirm dispensing amount.




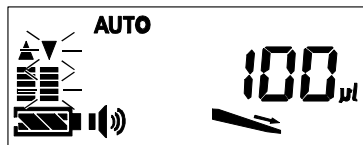
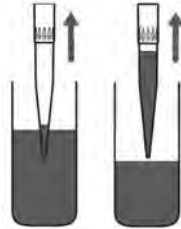
If you would like to change the dispensing amount, Press the  Enter key before aspirating.


When setting, operate from step 4 as described above.

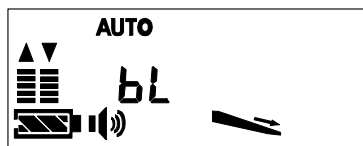
3) Operating the standard mode






[1] Put the tip end into the liquid to be dispensed and press the  Operation key on the pipette to aspirate the liquid.



[2] Put the tip end into the receiving vessel and press the  Operation key on the pipette to dispense the liquid.



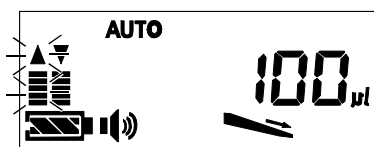[3] Press the  Operation key on the pipette to dispense the liquid remaining in the tip.
(When the blowout function is ON ().)




9-2 Multiple dispensing mode (MD)


1) Operation of the multiple dispensing mode


This is a suitable function to dispense the same liquid continuously, such as when dispensing a liquid on a microplate, etc. The operation consists of aspirating one time and dispensing several times. When carrying out pre-rinse for multiple dispensing or stopping multiple dispensing, use the total discharge function. (Refer to “9-8 Total discharge function”)

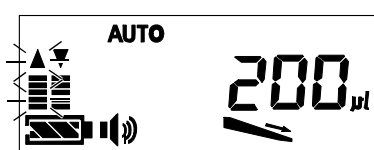2) Selecting the multiple dispensing mode





(The display is an example.)

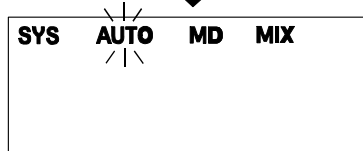
[1] Press the  Back key.



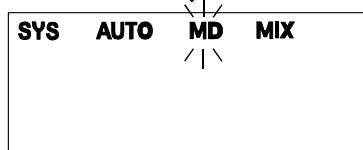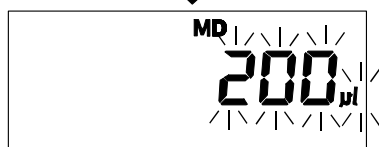
[2] Press the  Up or Down key to select "MD".




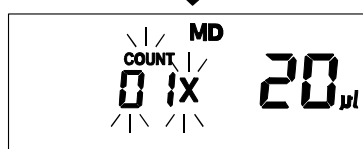
[3] Press the  Enter key to select multiple dispensing mode (MD).



[4] Press the  Up or Down key to set dispensing amount for one time.

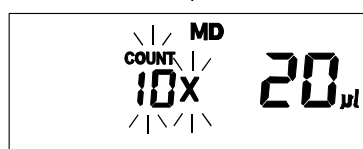



[5] Press the  Enter key to confirm dispensing amount for one time.





[6] Press the  Up or Down key to set dispensing count.

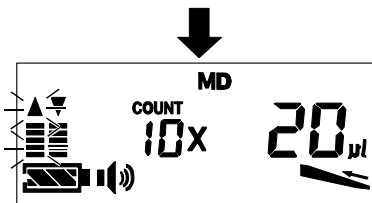* "Dispensing amount for one time x dispensing count" cannot be set if it exceeds the volume range.






[7] Press the  Enter key to confirm the dispensing count.


If you would like to change dispensing amount or dispensing count, press the  Enter key before starting aspiration.

When setting, operate from step 4 described above.



3) Operating the multiple dispensing mode

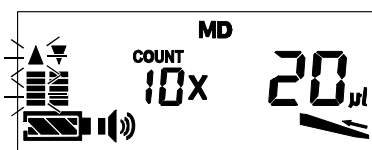
The following example is when dispensing 20 µ L x 10 times.






[1] Put the tip end into the liquid to be dispensed and press the  Operation key on the pipette to aspirate the liquid.

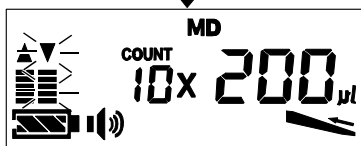
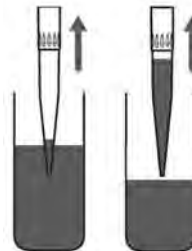
(The example shows 20 µ L x 10 times = approx. 200 µ L)



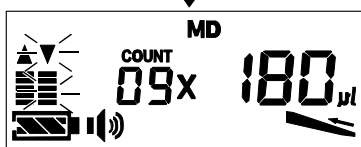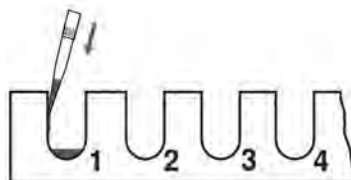[2] Put the tip end into a receiving vessel and press the  Operation key on the pipette to dispense the amount for one time.



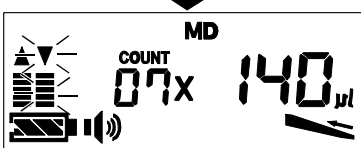
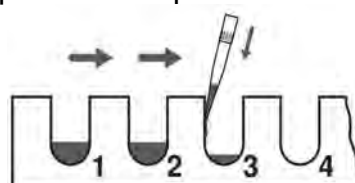
[3] In a similar manner, put the tip end into the next vessel and press the  Operation key to dispense the liquid.




[4] Repeat the operation described above with the dispensing count set.

When the set count of dispensing operations is completed a buzzer will sound twice.

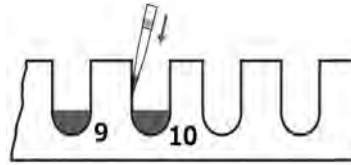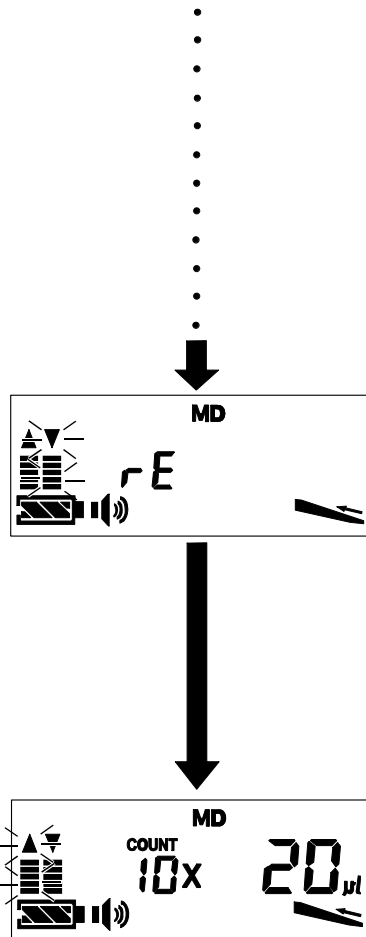(In the example, the buzzer sounds twice after 10 dispensing operations have been completed)



[5] Press the  Operation key on the pipette to discharge the remaining liquid.

(In the multiple dispensing mode, reverse operation occurs automatically.)

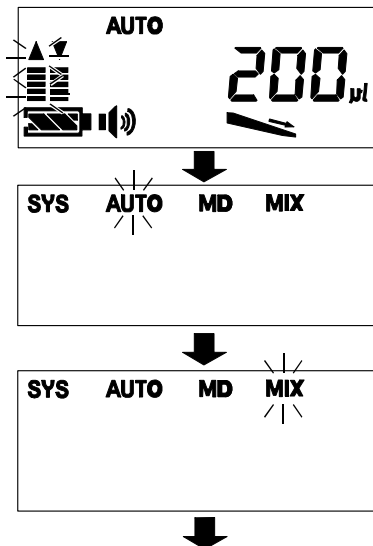


9-3 Mixing mode (MIX)

1) Operation of the mixing mode

This is a useful operation when uniformly mixing different types of liquids. In this method, aspirating and dispensing are repeated. This type of repetitive operation often results in fatigue, but with this pipette it is automatically carried out with the touch of one switch.


2) Selecting the mixing mode





(The display is example.)

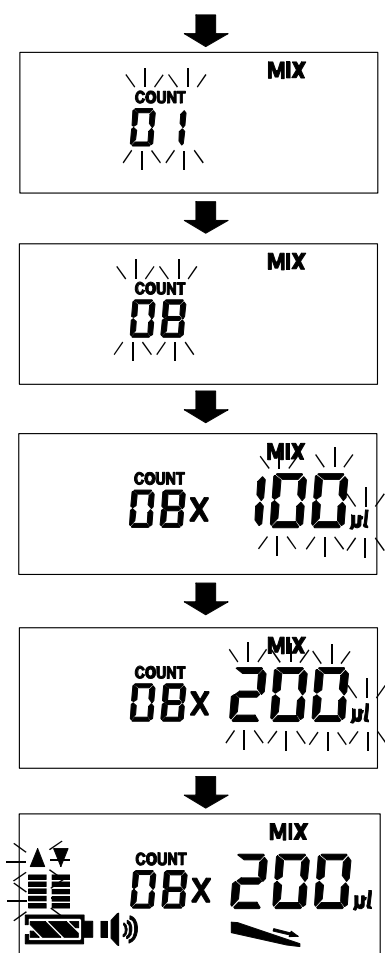[1] Press the  Back key.


[2] Press the  Up or Down key to select "MIX".





[3] Press the  Enter key to select mixing mode (MIX).

[4] Press the  Up or Down key to set mixing count (One set with aspirating and dispensing).
* Ten times maximum

[5] Press the  Enter key to confirm mixing count.

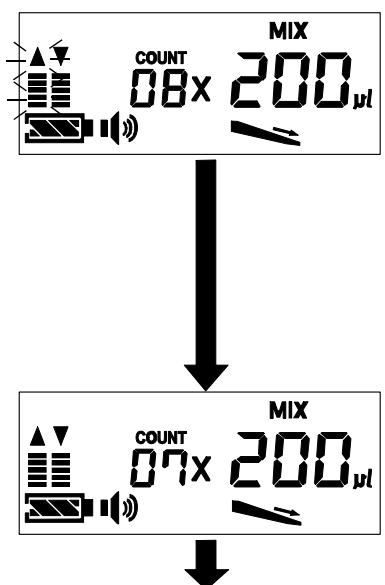
[6] Press the  Up or Down key to set mixing volume (Volume for aspirating when mixing).

[7] Press the  Enter key to confirm mixing volume.

If you would like to change the mixing count or amount, press the  Enter key before aspirating.


When setting, operate from step 4 described above.

3) Operating the mixing mode




[1] Insert the tip end in the liquid to be mixed.



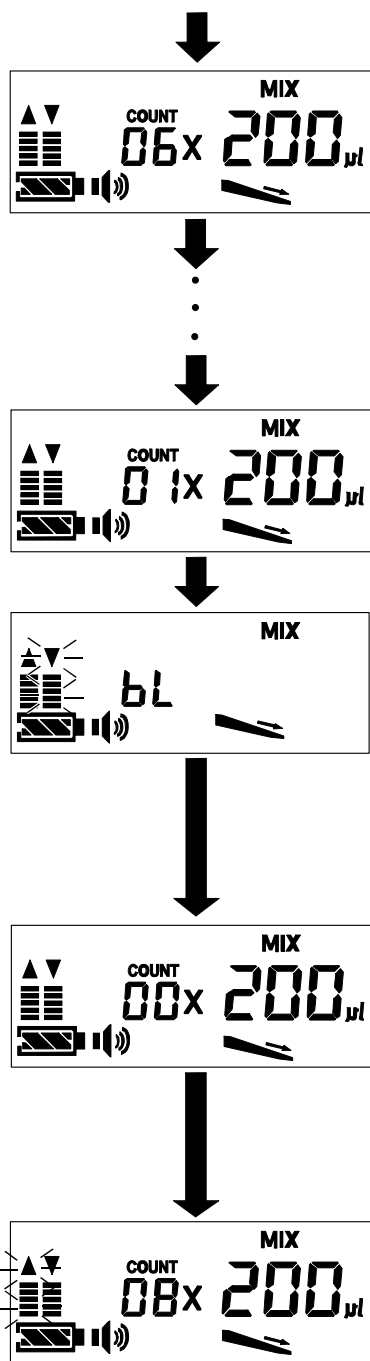
[2] Press the  Operation key on the pipette to aspirate the set mixing volume.





Aspirating of mixed

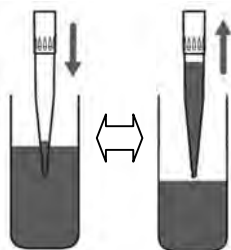

[3] The pipette executes the set count of aspirating-dispensing cycles by approx. 2/3 of the set mixing volume.



Dispensing to absorption
Repeat with mixing count set.

[4] Press the  Operation key on the pipette to dispense the liquid remaining in the tip.
(When the blowout function is ON () .)

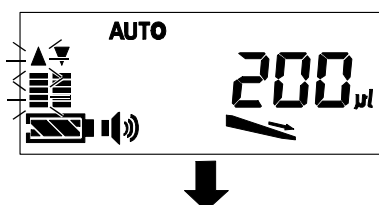


[5] One mixing operation is completed once the aspirated liquid is completely discharged.
When the operation is completed a buzzer sounds twice.

9-4 System setting mode (SYS)


The SYS mode is used to perform or read out pipetting operation settings that suit the purpose or the liquid to be handled.

2) Operating the system setting mode






(The display is example.)

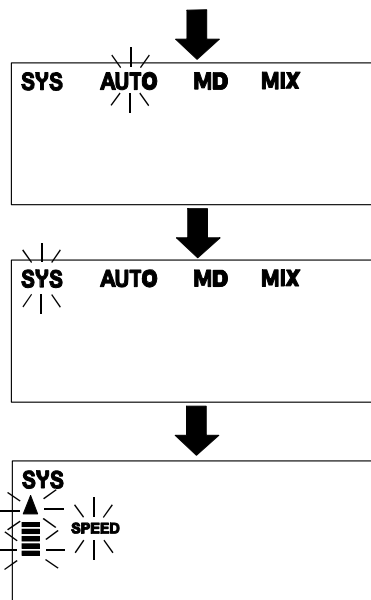[1] Press the  Back key.



















[2] Press the  Up or Down key to select “SYS”.

[3] Press the  Enter key to enter the system setting mode.

[4] By pressing the  Up or Down key, change the setting of each item, and press the  Enter key to confirm it.

2) Item of the system setting mode

Display of the each item and setting contents

Functions	Displays	Setting contents
Aspirating speed	 SPEED	High speed  Low speed 
Dispensing speed	 SPEED	High speed  Low speed 
Buzzer		On
		Off
Blowout		On
		Off
Reverse operation *3		Off
		On
Program memory	PROG  	Nine programs between 01 and 09 are available for the “   “. (Read out your preferred program from the programs previously set.)

*3 The reverse operation is only selectable when blowout setting is off. It cannot be selected when the pipette is in MIX mode, either.

9-5 Program setting mode

The nine programs can be saved in the program memory built into the pipette (PROG 01 to 09). By saving a frequently used mode or volume for operation beforehand, these setting can easily be read out from the next use. Select and set the mode or volume to be saved before saving the program setting.

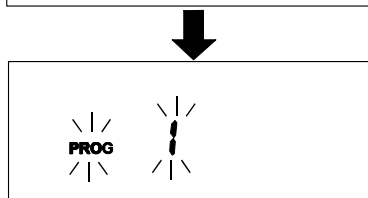
Saving the program setting


Set the pipette to your preferred settings.

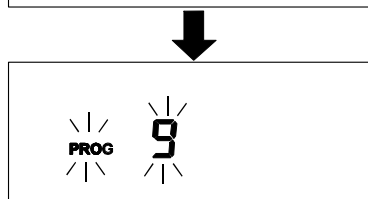
Example: When saving the AUTO mode, dispensing volume 200 μ L, buzzer ON and blowout ON




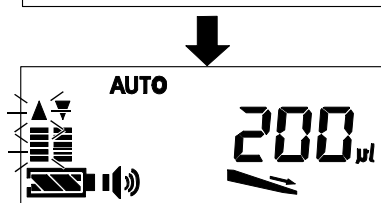
(The display is example.)




[1] Press and hold the  Down key to enter program setting saving mode.



[2] Press the  Up or Down key to select program number to be save from 1 to 9.



[3] Press the  Enter key to save the setting.


Once these have been saved a buzzer will sound once.

Reading out the program setting

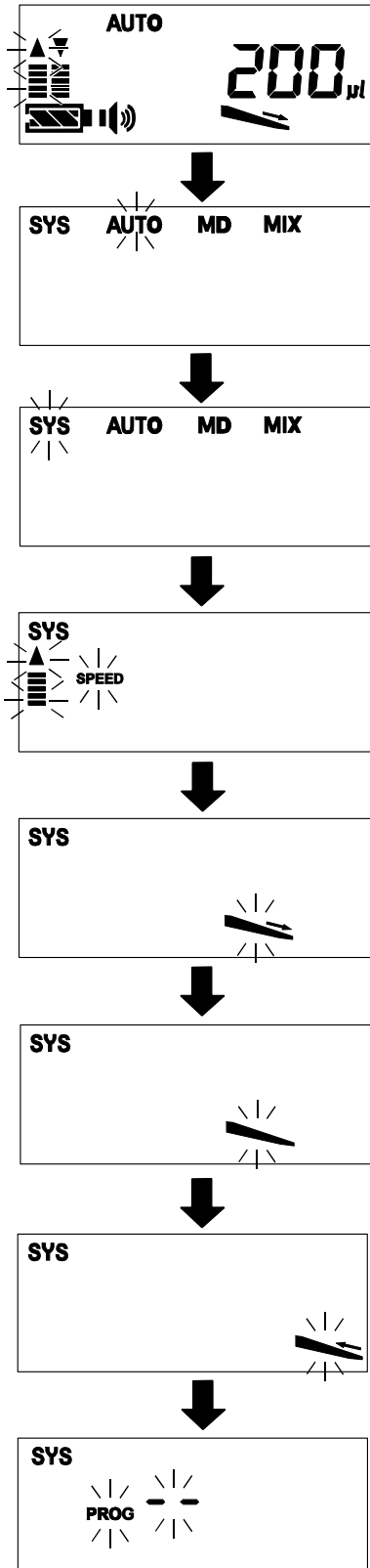
The set program can be read out at system setting mode (SYS).
(Refer to “9-4 system setting mode (SYS)” for details)

9-6 Reverse operation (Dispensing liquid that tends to remain in the tip)

When you would like to accurately dispense a highly viscous liquid that has a tendency to remain in the tip, we recommend using the reverse operation. By aspirating a large amount of the liquid beforehand, the reverse operation enables the correction of the amount of liquid remaining in the tip.

To enable reverse operation, set the setting of  at system setting mode (SYS).
(Refer to “9-4 System setting mode (SYS)”.)


1) Setting the reverse mode








(The display is example.)




[1] Press the  Back key.



[2] Press the  Up or Down Key to select "SYS".

[3] Press the  Enter key to enter the system setting mode.

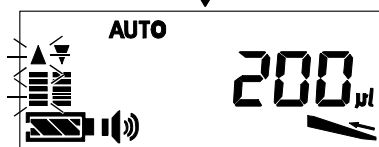
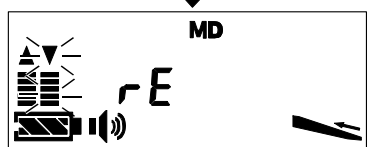
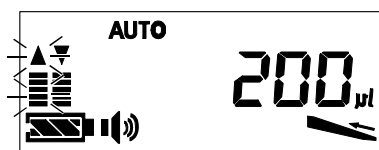
[4] Press the  Enter key several times until  is displayed.


[5] Turn the blowout operation mode OFF  using the  Up/Down key, then confirm this selection by pressing the  Enter key.

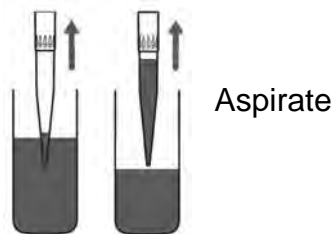
[6] Turn the reverse operation mode ON  using the  Up/Down key, then confirm this selection by pressing the  Enter key.


[7] End the system setting mode by pressing the  Enter key while **PROG**  is being displayed.

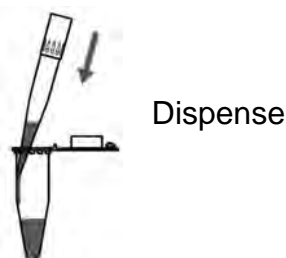
2) Operating the reverse mode




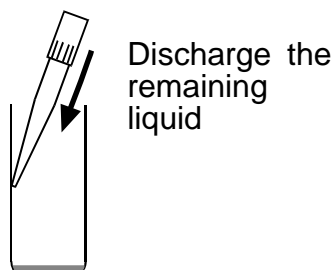
[1] Press the  Operation key on the pipette to aspirate the liquid.



[2] Press the  Operation key on the pipette to dispense the liquid.





[3] Press the  Operation key on the pipette to discharge the remaining liquid.




9-7 Blowout function


This is the function to forcibly dispensing the liquid remaining in the end of the tip by temporarily lowering the piston built in the pipette below the start position for aspiration after dispensing the liquid remaining in the tip.

By pressing the  Operation key when “bL” is shown on the display, carry out blowout.

* After carrying out blowout, the built-in piston remains in the blowout position while the  Operation key is being held down, and it returns to initial position when the finger was released from the Operation key. By releasing the Operation key after removing the tip end from the vessel, aspiration of the dispensed liquid in the tip again can be prevented.

9-8 Total discharge function

Pressing and holding the  Operation key expels all of the liquid remaining in the tip. This function is useful when, for example, you want to terminate the operation halfway through multiple dispensing.

Continuing to hold the  Operation key down after this turns the pipette power off.

9-9 “Dispensing correction function” for multiple dispensing

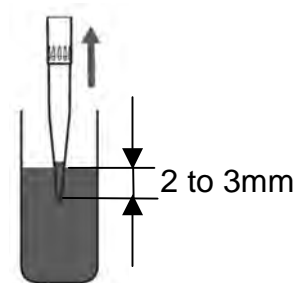
The electronic pipette aspirates and dispenses a liquid by moving the internal piston up and down using motor. Since movement of motor and piston reverses when the operation switches from aspirating to dispensing, an error in dispensing volume due to backlash will occur. To correct this error, the MPA series is equipped with the “dispensing correction function” (Patent pending) for multiple dispensing, which automatically discharges a small amount of sample before delivery. This ensures the piston is always set in the descending direction when dispensing starts, keeping the margin of error to a minimum.

10. PIPETTING FOR ACCURATE DISPENSING

- When performing aspiration, if the tip is immersed too deeply into the sample liquid, an amount larger than the selected dispensing volume may be delivered, as excess liquid attaches the outside of the tip. Ideally, for aspiration, the tip should be dipped into the liquid to a depth of 2 to 3 mm.


The pipette is designed to correctly perform aspiration when it is in the vertical position.

Therefore, hold the pipette as vertically as possible when aspirating.



- When replacing the tip, pre-rinse the tip with the necessary dispensing volume setting. The reverse operation is recommended for a sample liquid that tends to linger in the tip.

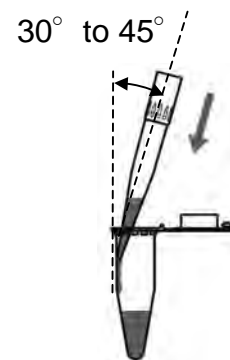- When dispensing the sample, touch the tip at a 30° to 40° angle to the wall of the receiving vessel.



11. CALIBRATING THE PIPETTE USING AN ELECTRONIC BALANCE

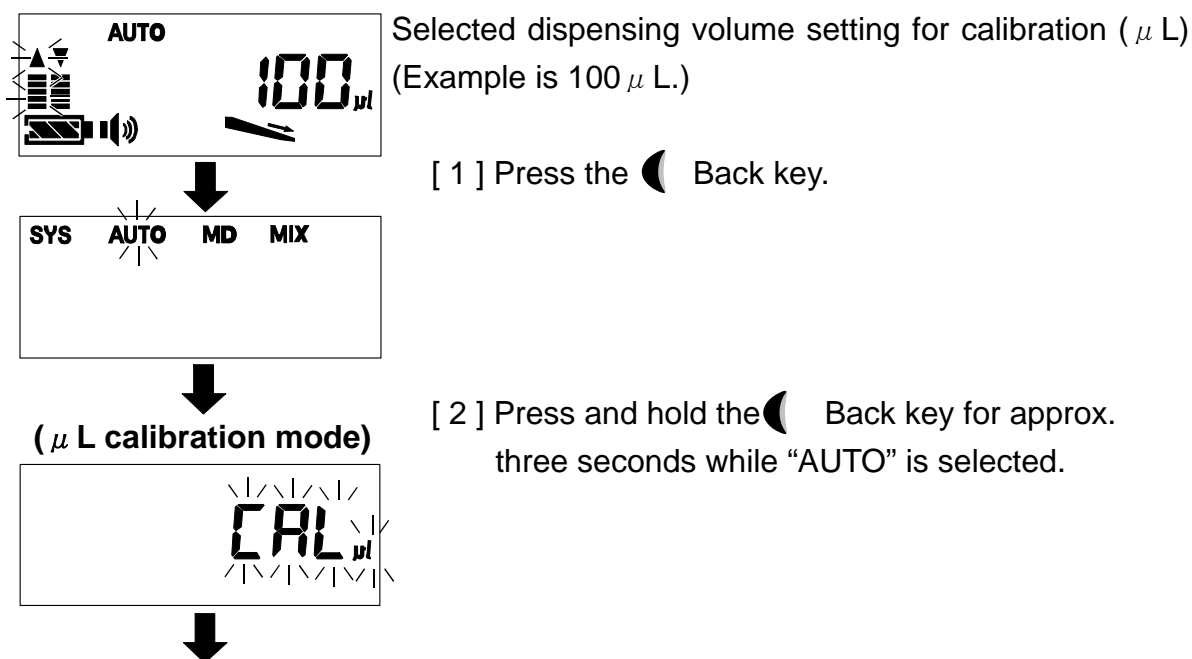
The MPA series provide user with a dispensing volume calibration function. Using this function, it is easy to correct (calibrate) errors due to differences in tips used, etc. When you need to always control the dispensed volumes in a precise manner, perform volume calibration as necessary when you change the dispensing volume setting. For verification of dispensed volumes necessary for calibration, A&D's pipette accuracy tester - AD-4212B-PT, FX-300i-PT, or combined use of BM series and BM-014 (Sold separately) - are useful.

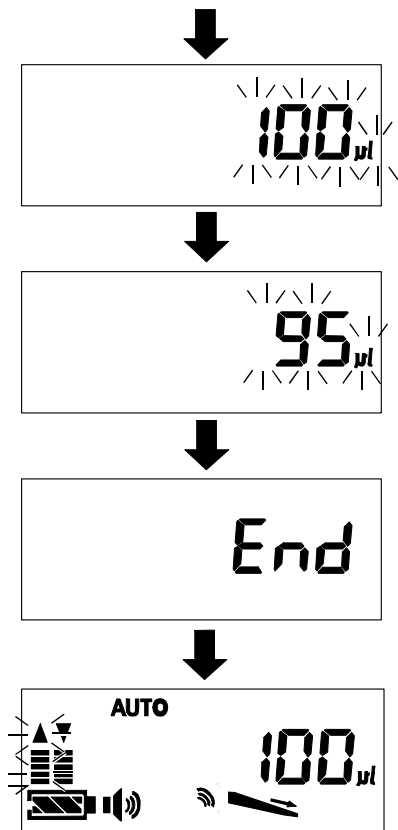
11-1 Volume calibration function (μ L calibration function)


This is a function to correct the dispensing volume of the MPA series. Using an A&D pipette accuracy tester or other appropriate device, measure the volume actually dispensed as opposed to the selected dispensing volume setting, and then enter the actual dispensed volume to the pipette to correct its dispensing volume. To calibrate the dispensing volume, complete the following procedure:


Calibrating the dispensing volume


1. Set the dispensing amount of the MPA series to the volume to which you would like to calibrate it. (The example is 100μ L)
2. Using an electronic balance, measure and record an actual dispensed volume as opposed to the selected dispensing amount setting. (The example is 95μ L)
3. Enter an actual dispensed volume to the pipette by the following procedure.






[3] Press the  Enter key to enter μ L calibration mode.

[4] Press the  Up or Down key to alter the value to the actual dispensed volume.

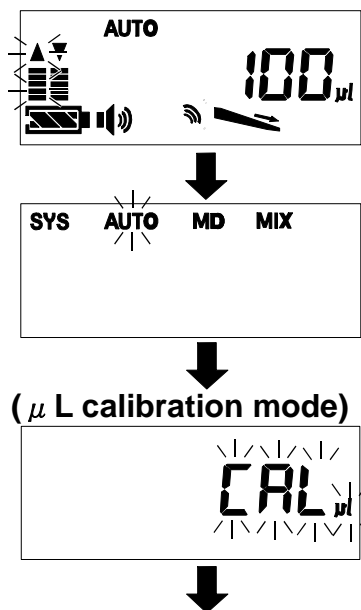
[5] Press the  Enter key to confirm actual dispensed volume.

[6] The notice mark  lights up to indicate that the volume calibration has been performed. (In this example, the dispensed amount is corrected and altered to 100 μ L.)

* After calibrating, a volume range that can be selected may be limited depending on available movement range of the piston.

11-2 Resetting the volume calibration


Go through the following procedure to restore the factory default settings for volume calibration:



(The display is example.)

[1] Press the  Back key.

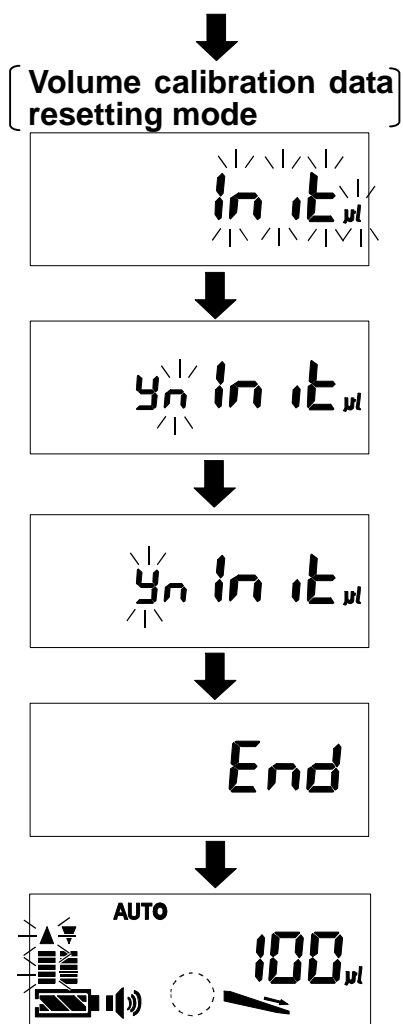[2] Press and hold the  Back key for approx. three seconds while "AUTO" is selected.



[3] Press the Up or Down key to select volume calibration data resetting mode.

[4] Press the Enter key to enter volume calibration data resetting mode.

[5] Press the Up or Down key to select **y** .
(**y** means YES and **n** means NO.)

[6] Press the Enter key to reset volume calibration data.

[7] The Notice mark is turned off to indicate that volume calibration data have reverted back to the factory default.

11-3 Dispensing in a unit of weight (in mg unit)

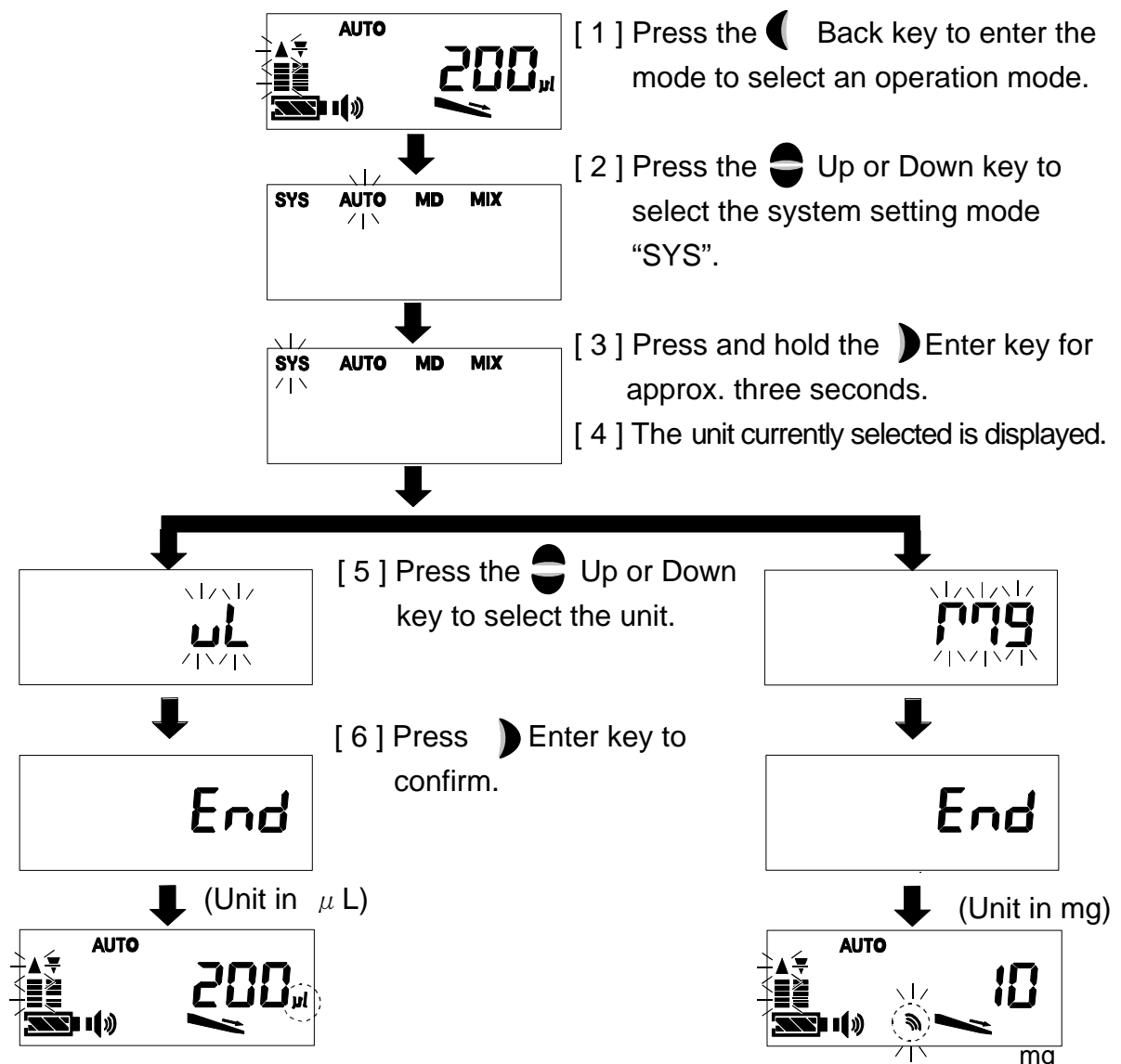
Dispensing of a liquid can be performed by weight (mg) instead of volume. (mg unit function)
This function is useful when you handle a liquid that needs to be managed by weight, such as a diluted solution of a solid or powder. Although the density of a liquid can vary depending on the sample type and concentration, by weighing the dispensed amount with an electronic balance and inputting the result into the pipette, it becomes possible to easily dispense the liquid in a unit of weight (mg).

Selecting the mg unit

The unit (volume: μL / weight: mg) for pipetting can be toggled by the following method.
When the mg unit is selected, the Notice mark blinks and the μL unit is turned off.

* When the unit of weight (mg unit) is selected, perform weight calibration by the dispensing amount to be used. (The weight calibration data reverts back to the factory default once the unit is switched to μL .)

Method for selecting the mg unit



With the μ L unit selected, the Notice mark is turned off while the μ L unit mark lights up.

With the mg unit selected, the Notice mark blinks while the μ L unit mark is turned off

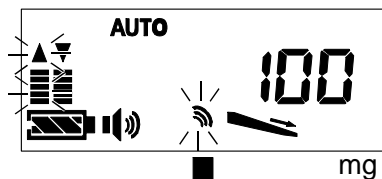
* The calibration data reverts back to the factory default once the weight unit is switched.


11-4 Weight calibration function (mg calibration function)


The density of a liquid varies depending on the type and concentration of the material. Make sure to perform mg calibration when you dispense a different sample or use the mg unit for the first time. Further, when you need to always control dispensed amounts in a precise manner, perform mg calibration when you change the dispensing amount setting as well.

Method for mg calibration

1. Select the mg unit beforehand. (Refer to “11-3 Dispensing in a unit of weight (in mg unit)”)
2. Set the dispensing amount of the MPA series to the weight to which you would like to calibrate it. (The example is 100mg)
3. Using an electronic balance, measure and record an actually dispensed weight as opposed to the selected dispensing amount setting. (The example is 95mg)
4. Enter an actual dispensed weight to the pipette by the following procedure.



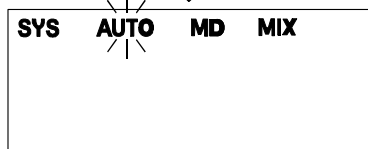
[1] Press the  Back key to enter the mode to select an operation mode.




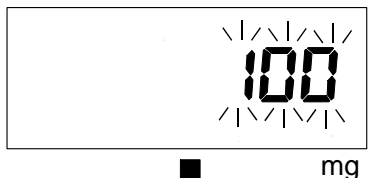
[2] Press and hold the  Back key for approx. three seconds while “AUTO” is selected.


(mg calibration mode)

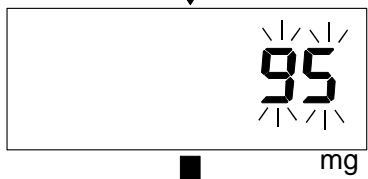



[3] Press the  Enter key to enter mg calibration mode.


In mg calibration mode, the Notice mark  illuminates while μ L mark is turned off




[4] Press the  Up or Down key to alter the value to the actual dispensed weight. (The example is 95mg)



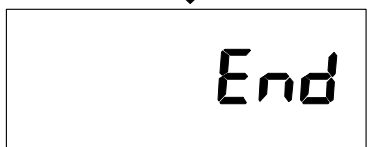[5] Press the  Enter key to confirm the actual dispensed weight.



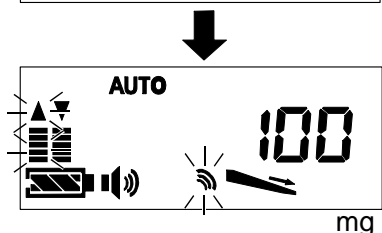(In this example, the dispensed amount is corrected and altered to 100mg.)




12. ADJUSTING HEIGHT OF THE TIP EJECTOR

A height of the tip ejector can be adjusted so that it can match the conditions of how the tip used was connected.


Use a small minus screwdriver.

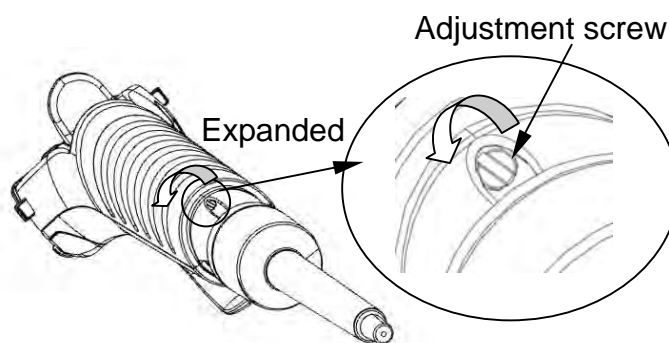By turning the adjustment screw in a counter-clockwise direction, the tip ejector can be lowered toward the tip side.



13. STORAGE AND MAINTENANCE

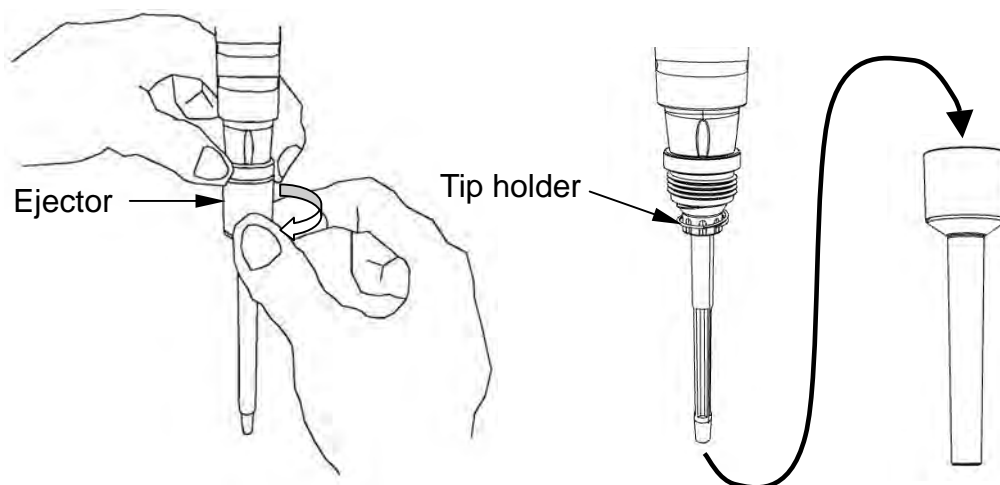When cleaning the outside of the pipette or the tip holder, remove dirt by using a cloth dampened with 60% isopropyl, 70% ethanol or ph-neutral detergent, and wipe it using a lint-free and dry soft cloth. Also, because the ejector can easily be removed, the tip holder can easily be cleaned.



Regularly clean the tip holder about once every half a year to a year, depending on the frequency used. Regularly confirm the pipette accuracy.

Maintenance

When inspecting the pipette or managing pipette precision, the following device can be used. (Refer to “17. LIST OF ITEMS SOLD SEPARATELY (DISPOSABLE ITEMS) ”)

- Leak tester (Sold separately: AD-1690):

Leakage within the pipette can be easily checked.

- Pipette accuracy tester (Sold separately: BM series, AD-4212B-PT, FX-300i-PT):

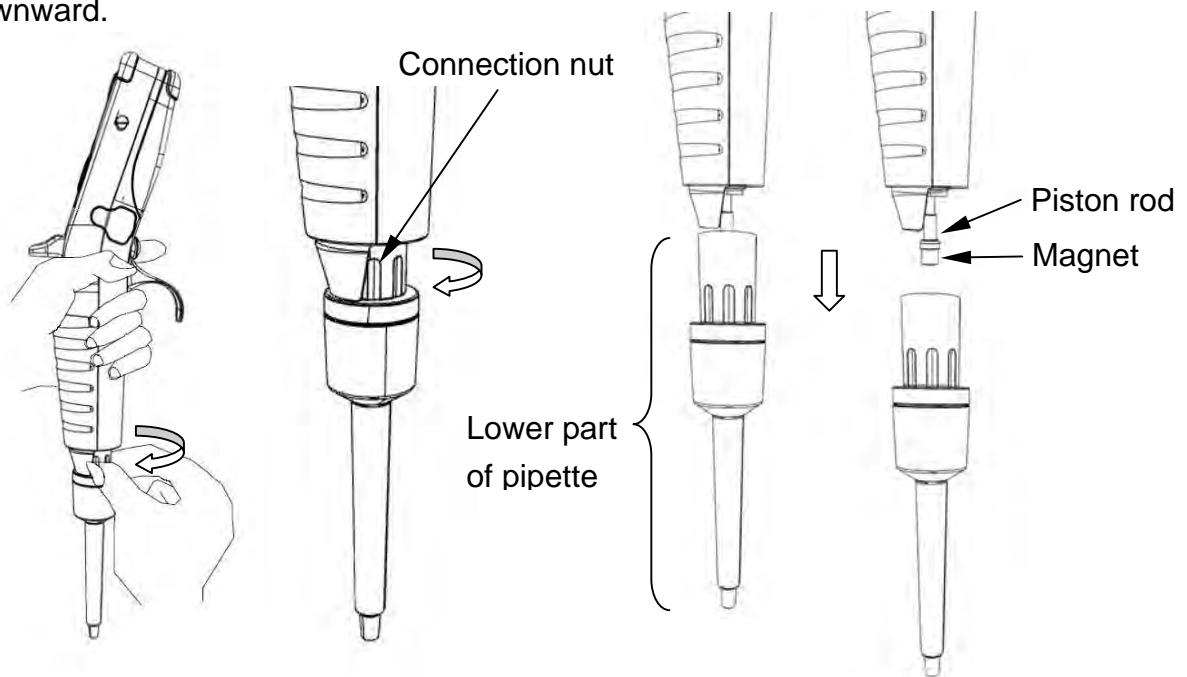
By using an electronic balance and specified software “WinCT-pipette”, the actual dispensed amount of the pipette can be confirmed.

13-1 Autoclave

By removing the lower part of pipette from the pipette, it can be placed in an autoclave.
Condition for autoclave: Lower part of pipette 121°C 2 atm 20 minutes

Removing the lower part of pipette

- The piston and piston rod are connected using a magnet. Therefore, by removing the connection nut completely, the lower part of pipette can be removed by pulling downward.



- The pipette has a strong magnet. Therefore, do not allow any magnetic material to be close to the magnet, or do not allow an instrument that is affected by magnetic force to be close to the magnet.
- Do not autoclave the lower part of the pipette beyond the stated specifications as this will result in damage to the pipette.
- After performing autoclave, install the lower part of the pipette to the pipette after drying it completely. Keep the lower part of the pipette in a room for at least two hours before installing.

14. TROUBLE SHOOTING

The electronic pipette has a self-check function. When there is an error in these results, the pipette displays a following error message “Err XX”. When an error cannot be released even after applying a following remedy to the pipette, please apply for a repair. (Refer to “15. WHEN REQUESTING REPAIR”.)

Circumstances	Considerations for cause	Remedies
“Err 01”	The automatic position resetting is not working.	If the connection nut (Refer to “7-5 Parts names”) is loose, retighten it. Press the Reset key to reset the pipette.
“Err 02”	Step motor error.	Press the Reset key to reset the pipette.
The sample liquid lingers inside the pipette tip.	Compatibility in characteristics between the sample liquid and pipette tip.	Select the reverse operation. Use a tip that matches the sample liquid.
The sample liquid leaks from the pipette tip, or the aspirated volume is small.	A wet tip was used.	Use a new tip.
	The tip is improperly attached.	Attach the tip properly.
	Quality error of the pipette tip.	Use a high quality tip.
	Sealing problem in the piston area.	Replace the lower part of the pipette.
	Wear on the tip holder.	Replace the tip holder.
Aspirating error	Weak battery.	Recharge the battery.
	The tip ejector is improperly attached.	Disassemble and then reassemble the lower part of the pipette.
	The tip holder is clogged.	Remove the clog.
	The piston movement is not smooth.	Clean the piston. Replace the lower part of the pipette
The power cannot be turned on.	Contact failure of the battery.	Set the battery again.
	Dirt on the battery terminal.	Replace the battery with a new one.
	Battery discharge.	Recharge the battery.
	Deterioration of the battery.	Replace the battery with a new one.

15. WHEN REQUESTING REPAIRS

The pipette requires repair if an error occurs and cannot be corrected by following the troubleshooting methods provided in this manual. In this case, please contact your local A&D representative.

When requesting repairs, it is essential that you confirm the pipette is free of contamination by a harmful material. Please photocopy the "Attestation of contamination removal" that can be found on the last page of this manual, fill in the required items, and attach it to the pipette you are going to send.

16. SPECIFICATION

		MPA-10		MPA-20		MPA-200		MPA-1200	
Volume range		0.5 to 10.0 μ L		2.0 to 20.0 μ L		10 to 200 μ L		100 to 1200 μ L	
Performance *1 *2	Volume	1.0 μ L	10.0 μ L	2.0 μ L	20.0 μ L	10 μ L	200 μ L	100 μ L	1200 μ L
	Accuracy	\pm 4.0%	\pm 1.0%	\pm 4.0%	\pm 1.0%	\pm 2.5%	\pm 0.6%	\pm 2.5%	\pm 0.5%
	Repeatability (CV)	2.5%	0.4%	2.5%	0.4%	1.0%	0.15%	0.6%	0.15%
Operation mode		Standard mode (AUTO), Multiple dispensing mode (MD), Mixing mode (MIX), System setting mode (SYS)							
Program memory		9 programs							
Aspirating and dispensing speed		5 speed							
Pipette driving method		High performance stepping motor							
Energy saving setting		Automatically power turning off after ten minutes							
Maximum number of dispensing (When recharging fully)		Approx. 1,800 times *2							
Charging time		Approx. 5 hours at 100% charging							
AC adapter *3		- Input: AC100-240V 50/60Hz Power plug: Selectable - Output: DC5V / 1A							
Autoclave processing		Available only for the lower part of the pipette (121°C, 2 atm, 20 minutes)							
Use environment temperature		15 to 30°C							
Use environment humidity		RH 85% or less							
Battery		Lithium-ion battery 3.7V / 920mAh							
Total length		Approx. 280mm							
Weight (Battery is included.)		Approx. 150g				Approx. 160g		Approx. 170g	

*1 Performance with the MPA-10 is ensured with 1.0 μ L or more.

*2 When setting to standard mode, maximum aspirating and dispensing speed.

*3 For recharging. The pipette can be used even when recharging.


Note: The specification may be changed without notice.

17. LIST OF ITEMS SOLD SEPARATELY (DISPOSABLE ITEMS)


17-1 Stands and hanger

Number	Name	Applicable pipette			
		MPA-10	MPA-20	MPA-200	MPA-1200
AX-ST-CHG	Charging stands	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
AX-HA-CHG	Charging hanger	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
AX-ST-ACR	Acrylic stands	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
AX-ST-SUS	Stainless steel stands	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
AX-HA-STD	Hanger	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>

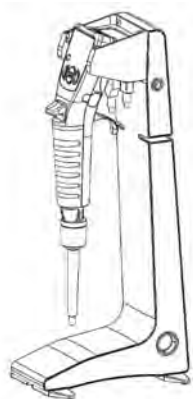
AX-ST-CHG Charging stands




AX-ST-CHG Charging hanger

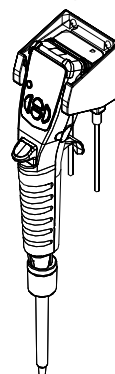


AX-ST-ACR Acrylic stands

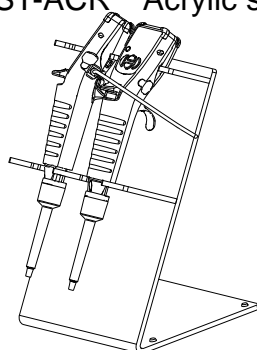

AX-ST-SUS Stainless steel stands

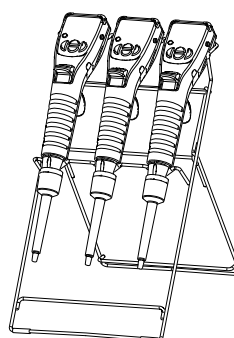

AX-HA-STD Hanger

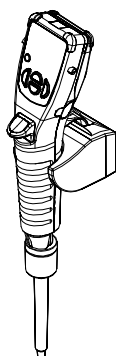

17-2 Tips and accessories

Number	Name	Applicable pipette			
		MPA-10	MPA-20	MPA-200	MPA-1200
AX-BOX-200A	Tip box with locking parts *4 *5 (10/20/200 μ L)	○	○	○	
AX-BOX-1200A	Tip box with locking parts *4 *5 (1200 μ L)				○
AX-BOX-200B	Tip box *4 *6 (10/20/200 μ L) \times 10 set	○	○	○	
AX-BOX-1200B	Tip box *4 *6 (1200 μ L) \times 10 set				○
AX-CART-10/20	Tip cartridge *7 *8 (10/20 μ L) \times 10 set A&D 10/20 μ L Standard tip	○	○		
AX-CART-200	Tip cartridge *7 *8 (200 μ L) \times 10 set A&D 200 μ L Standard tip			○	
AX-CART-1200	Tip cartridge *7 *8 (1200 μ L) \times 10 set A&D 1200 μ L Standard tip				○

*4: The tip is not included with the tip box.

*5: Material: Upper case, lower case and locking parts...PC

Base...TPE

Label...PET

(Available for autoclave performing)

* When performing autoclave to the pipette,
remove the base from the tip box.

*6: Material: Upper case and lower case ...PC

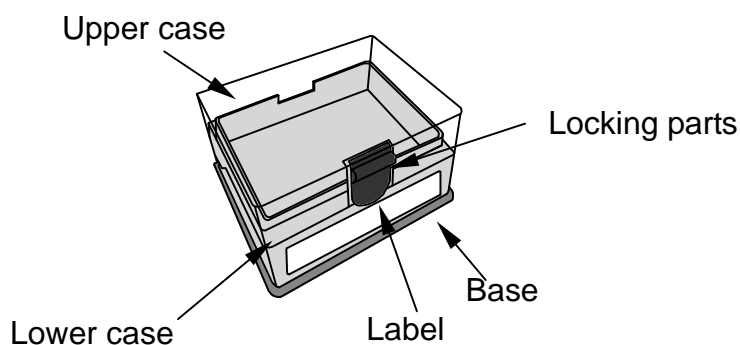
(Available for autoclave performing)

*7: 1 set = 96 tips

*8: Material: Tip, cartridge...PP


Example)

○AX-BOX-1200A Tip box with locking parts




○AX-BOX-1200B Tip box

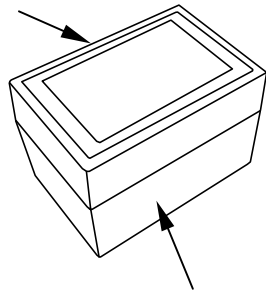Upper case




Lower case

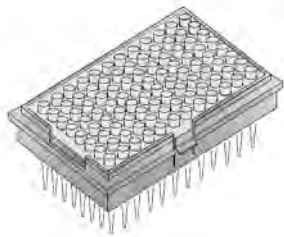○AX-CART-1200 Tip cartridge (1200 μ L)



Consultation)

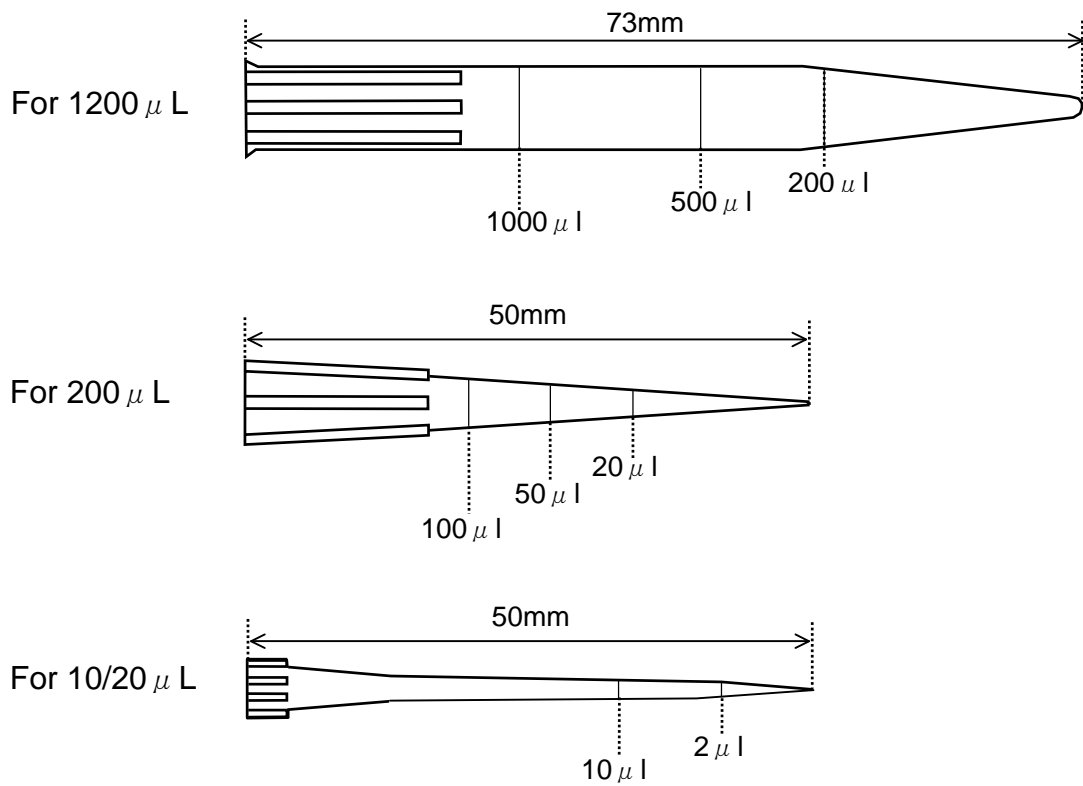The tip is marked with lines to act as guidelines for the aspiration amount.



17-3 Disposable items (User replaceable)



Number	Name	Applicable pipette			
		MPA-10	MPA-20	MPA-200	MPA-1200
AX-LOW-10	Lower part (10 μ L)	○			
AX-LOW-20	Lower part (20 μ L)		○		
AX-LOW-200	Lower part (200 μ L)			○	
AX-LOW-1200	Lower part (1200 μ L)				○
AX-BAT-MPA	Battery	○	○	○	○
AX-TB265	AC adapter (Provided as standard)	○	○	○	○



Lower part



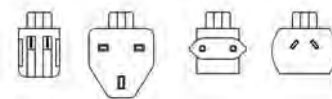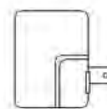Battery




AC adapter

17-4 Others

- Leak tester AD-1690

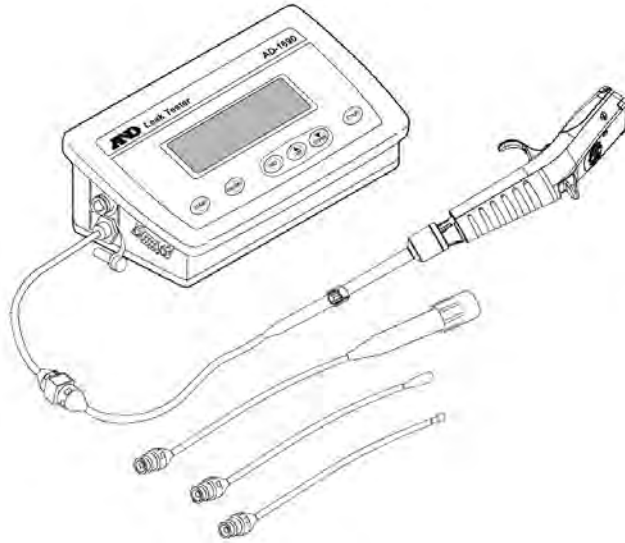
Leakage within the pipette can be easily checked.




AD-1690

- Pipette accuracy tester

{	BM-20/22(BM-014 attached)	10 μ L/20 μ L
	AD4212B-PT	10 μ L/20 μ L/200 μ L
	BM-252(BM-014 attached)	10 μ L/20 μ L/200 μ L/1200 μ L
	FX-300i-PT	1200 μ L

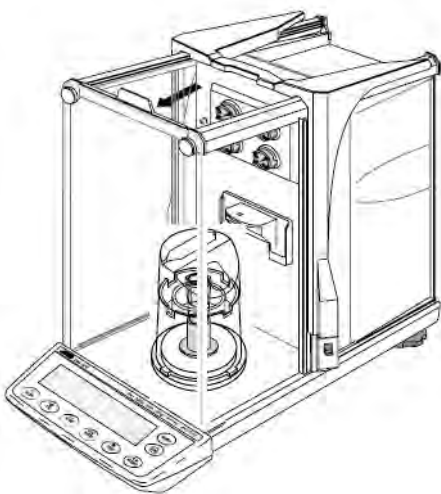

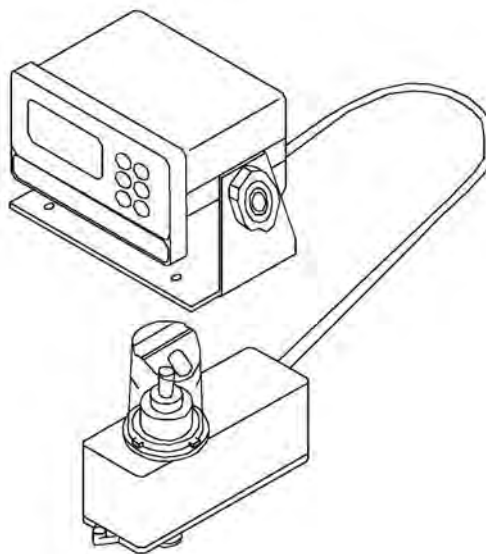BM-20/22 (BM-014 attached)
BM-252 (BM-014 attached)



AD-4212B-PT



FX-300i-PT

Attestation of contamination removal

Please fill in the following items when you send a pipette for repair.

Model name: _____

Serial number S/N: _____

I attest to the fact that this pipette is free of contamination by any substances that could pose a health threat to humans, such as Infectious bacteria or viruses, radioactive substances with associated risks of exposure, toxic substances, etc.

Signature: _____

Date: _____

Company name

(Facility name): _____

Section name: _____

Address: _____



A&D Company, Limited

3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013, JAPAN
Telephone: [81] (3) 5391-6132 Fax: [81] (3) 5391-6148

A&D ENGINEERING, INC.

1756 Automation Parkway, San Jose, California 95131, U.S.A.
Telephone: [1] (408) 263-5333 Fax: [1] (408)263-0119

A&D INSTRUMENTS LIMITED

Unit 24/26 Blacklands Way, Abingdon Business Park, Abingdon, Oxfordshire OX14 1DY United Kingdom
Telephone: [44] (1235) 550420 Fax: [44] (1235) 550485

A&D AUSTRALASIA PTY LTD

32 Dew Street, Thebarton, South Australia 5031, AUSTRALIA
Telephone: [61] (8) 8301-8100 Fax: [61] (8) 8352-7409

A&D KOREA Limited

한국에이.엔.디(주)
서울특별시 영등포구 국제금융로6길33 (여의도동) 맨하탄빌딩 817 우편 번호 150-749
(817, Manhattan Bldg., 33. Gukjegeumyung-ro 6-gil, Yeongdeungpo-gu, Seoul, 150-749 Korea)
전화: [82] (2) 780-4101 팩스: [82] (2) 782-4280

ООО A&D RUS

ООО "ЭЙ энд ДИ РУС"
121357, Российская Федерация, г.Москва, ул. Вереysкая, дом 17
(Business-Center "Vereyskaya Plaza-2" 121357, Russian Federation, Moscow, Vereyskaya Street 17)
тел.: [7] (495) 937-33-44 факс: [7] (495) 937-55-66

A&D INSTRUMENTS INDIA PRIVATE LIMITED

ऐ&डी इन्सट्रुमेन्ट्स इण्डिया प्रा० लिमिटेड
509, उद्योग विहार , फेस -5, गुडगांव - 122016, हरियाणा , भारत
(509, Udyog Vihar, Phase-V, Gurgaon - 122 016, Haryana, India)
फोन : 91-124-4715555 फैक्स : 91-124-4715599